

EXPERTS IN PATIENT HANDLING SOLUTIONS




# NON-SLIP CHAIR PAD

This guide covers all products with the prefix NSP in the product code.

A low-profile breathable non-slip product which prevents slipping down in a chair or wheelchair. It does not help with repositioning. For repositioning products please use our **Glide and Lock Range**.


## USER INSTRUCTIONS

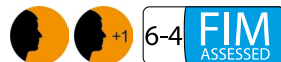

1

Place the **Non-slip Chair Pad** centrally on the chair/wheelchair with the breathable surface (BLUE) on the top and the non-slip base (BLACK) against the chair surface.

Check that the brakes are secure on the chair before commencing a movement. Wherever possible ask the person to help move themselves.

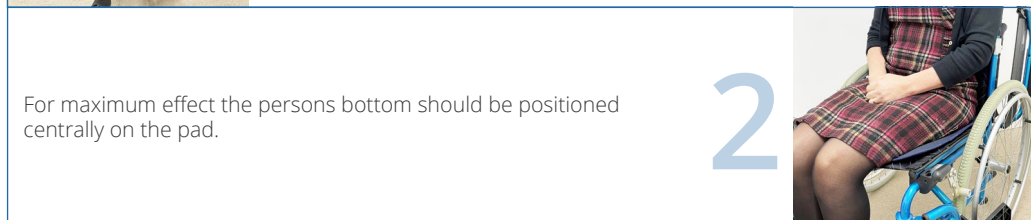


2

For maximum effect the persons bottom should be positioned centrally on the pad.



3

For some people it may also help to place a **Anti-slip Foot Pad** on the floor under their feet for added grip and/or patient safety.

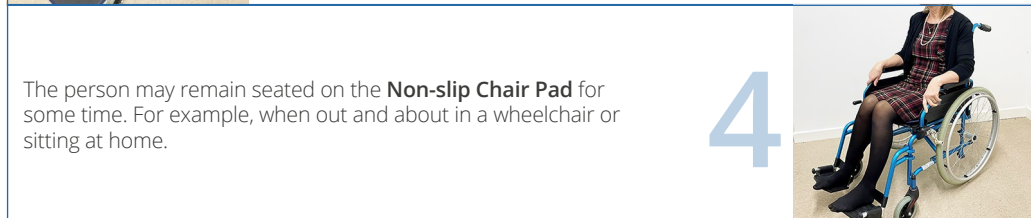

4

The person may remain seated on the **Non-slip Chair Pad** for some time. For example, when out and about in a wheelchair or sitting at home.

BREATHABLE TOP  
SURFACE (BLUE)



SECURE NON-SLIP  
BASE (BLACK)



## FABRIC CONTENT

Base fabric (Black): 100% Polyester with 100% PVC coating  
Top Fabric (Blue): 100% Breathable Polyester  
Filling: 100% PU Foam

## STORAGE AND HANDLING

- Store in a clean dry environment that is not subject to excessive heat.
- Keep away from sunlight.
- Store flat and avoid placing heavy objects on top of the product.



Keep dry



Keep away from  
sunlight

## WARNINGS AND PRECAUTIONS


- Do not leave on the floor when not being used.
- Do not use on a wet floor .
- Inspect before use for signs of wear or damage.
- Always test product grip on your chosen surface (floor, bed, driveway) before using.



Caution



Do not use if  
package is damaged

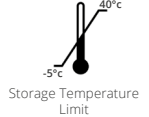

Storage Temperature  
Limit

## CLEANING INSTRUCTIONS

- Wipeclean with natural water and detergent or wipes that don't contain alcohol, solvents, bleaching or abrasive agents.
- Cleaning materials used should be patient safe and biodegradable.



Wipeclean only

## DISPOSAL INSTRUCTIONS

- Non-recyclable.
- Dispose of as clinical waste if it has been used by an infectious person.
- Otherwise dispose of through normal waste management.



Patient Information  
website



Medical Device



### MANUFACTURER

Hospital Direct (Marketing) Ltd, Units 3 & 4, The Green,  
Clun, Shropshire SY7 8LG, United Kingdom

**EU REPRESENTATIVE:** Med-Di-Dia Limited, Galway Technology Centre,  
Mervue Business Park, Galway H91 D932, Ireland



[www.hospitaldirect.co.uk](http://www.hospitaldirect.co.uk)

## HOSPITAL DIRECT (HD) ADDITIONAL SYMBOLS


At HD we understand that not everyone is familiar with all equipment, so to help therapists and healthcare professionals in assessing a product's suitability for a client, we have included two assessment tools, Functional Independence Measure (FIM) and Easy Guide Symbols. These will enable the healthcare professional to assess their client and decide whether or not a product will help the person based on their level of dependence and functional ability.

### FUNCTIONAL INDEPENDENCE MEASURE (FIM) ASSESSMENT

Our FIM assessment guide for each product will help you decide the suitability of this product for the person's ability and need. Based on the standard criteria from Level 6 where the person can use the product unaided and unsupervised to Level 1 where all the assistance is provided by the carer and the client can do nothing, this guide is easy, quick and simple to use to check suitability against ability and circumstance.

### EASY GUIDE SYMBOLS

We use four symbols to indicate where a product is suitable and safe for the person to use unassisted, with a single carer or with multiple carers. These symbols indicate the minimum recommendation.

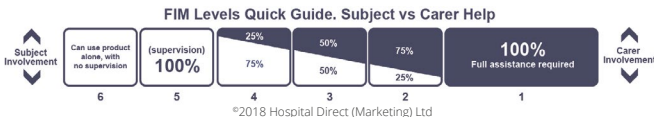

Person  
Unassisted



Single  
Carer



Two  
Carers



Four  
Carers

HD has a policy of continuous development and, as such, reserves the right to alter specifications (including measurements, materials and colours) without prior notice.